## Corning<sup>®</sup> PuraMatrix<sup>™</sup> Peptide Hydrogel

Catalog No. 354250

**Guidelines for Use** 

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences



### **TABLE OF CONTENTS**

INTENDED USE	2
MATERIALS PROVIDED	3
SUPPLEMENTAL/OPTIONAL MATERIALS (NOT SUPPLIED)	3
PRECAUTIONS	3
GENERAL USE	4
<ul> <li>PROCEDURES FOR USE</li></ul>	4 5 6 8 8 12 13 14
STORAGE	17
TECHNICAL SUPPORT	16
CUSTOMER SERVICE	16

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences



Page

For Research Use Only. Not for use in diagnostic or therapeutic procedures. For a listing of trademarks, visit www.corning.com/lifesciences/trademarks

© 2013 Corning Incorporated

### INTENDED USE

Corning<sup>®</sup> PuraMatrix<sup>™</sup> Peptide Hydrogel (Corning PuraMatrix) is a synthetic matrix that is used to create defined three dimensional (3D) microenvironments for a variety of cell culture experiments. To achieve optimal cell growth and differentiation, it is necessary to determine the appropriate mixture of this material and bioactive molecules (e.g., growth factors, extracellular matrix (ECM) proteins, and/or other molecules). Corning PuraMatrix consists of standard amino acids (1% w/v) and 99% water. Under physiological conditions, the peptide component of Corning PuraMatrix self-assembles into a 3D hydrogel that exhibits a nanometer scale fibrous structure with an average pore size of 50-200 nm. The hydrogel is readily formed in a culture dish, multiwell plate, or cell culture insert.

Corning PuraMatrix provides the benefits of 3D cell culture while allowing more straightforward 2D plating. Many cell types form appropriate 3D morphologies when plated on the surface of Corning PuraMatrix, which simplifies the work required to develop 3D *in vitro* model systems compared with cell encapsulation or sandwich culture systems. Research has shown that cells forming 3D structures often resemble their *in vivo* counterparts more closely when compared to cells grown in 2D. Corning PuraMatrix allows defined microenvironments to better recapitulate the *in vivo* milieu, which should improve the predictive quality of results obtained from assays conducted in a high-throughput context.

The hydrogel has been shown to support the differentiation of hepatocyte progenitor cells<sup>1</sup>, rat pheochromocytoma cells (PC12)<sup>2</sup> and hippocampal neurons<sup>3</sup>, as well as enabling endothelial cell tubulogenesis<sup>4</sup>. Studies have also demonstrated that Corning PuraMatrix supports the attachment of a variety of primary (e.g., neuronal, fibroblast, and keratinocyte) and transformed (e.g., MG-63, SH-SY5Y, HEK293, and NIH3T3) cell types<sup>5,6</sup>. Other potential applications include stem cell proliferation and differentiation, tumor cell migration, and *in vivo* analyses of tissue regeneration. Corning PuraMatrix is biocompatible, resorbable, and devoid of animal-derived material and pathogens. For *in vivo* studies in animals, the soluble material can be injected and will subsequently form a 3D hydrogel upon contact with the physiological environment.

### References

- 1. Semino, C.E., et al., *Differentiation* 71:262 (2003).
- 2. Holmes, T.C., et al., PNAS USA 97:6728 (2000).
- 3. Semino, C.E., et al., *Tissue Engineering* **10**:643 (2004).
- 4. Narmoneva, D., et al., *Biomaterials* 26:4837 (2005).
- 5. Zhang, S., et al., *Biomaterials* **16**:1385 (1995).
- 6. Thonhoff, J.R., et al., *Brain Research* **1187**:42 (2008).

For additional information, visit our website at: <a href="http://www.corning.com/lifesciences">www.corning.com/lifesciences</a>

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences



For Research Use Only. Not for use in diagnostic or therapeutic procedures. For a listing of trademarks, visit www.corning.com/lifesciences/trademarks © 2013 Corning Incorporated

2

### MATERIALS PROVIDED

• Corning<sup>®</sup> PuraMatrix<sup>™</sup> Peptide Hydrogel (Cat. No. 354250) packaged in one vial containing 1% solution (w/v) of purified synthetic peptide (5 ml per vial)

### SUPPLEMENTAL/OPTIONAL MATERIALS (NOT SUPPLIED)

- Tissue culture media for cells of interest
- 20% sucrose in distilled water (dH<sub>2</sub>O), sterile filtered (prepare 10% sucrose stock by diluting with sterile water)
- ECMs, growth factors, and cytokines
- Falcon<sup>®</sup> Cell Culture Inserts and Falcon Cell Culture Plates
- For information about available products from Corning Life Sciences, visit our website at: www.corning.com/lifesciences

### Hints for Optimal Results with Corning PuraMatrix Peptide Hydrogel

### PRECAUTIONS

- **GELATION:** Since the gelation of this material is initiated by salt concentrations ≥1 mM, this process may be promoted by salt-containing buffers and/or media. Therefore, do not combine Corning PuraMatrix with salt-containing buffers and/or media until gelation is desired. Once the gel is formed, any subsequent mixing of the sample will disrupt the structure of the hydrogel.
- **VISCOSITY:** Decrease the viscosity of the stock solution (1% w/v) by vortexing or sonication for 30 minutes in a bath sonicator every time the stock is used. If air bubbles are present, centrifuge an aliquot of the stock in a 1.5 ml microtube at full speed.
- **pH and OSMOLARITY:** Since Corning PuraMatrix exhibits a pH of 2-2.5, work quickly to minimize the amount of time that cells are in contact with the material prior to the addition of media when performing cell encapsulation experiments. Addition of isoosmotic sucrose to the cell suspension as well as to the Corning PuraMatrix, will help protect the cells until the pH is normalized by equilibration with tissue culture media. For details, see **Parts D and E (Cell Encapsulation Protocols)**.
- HANDLING:
  - Since Corning PuraMatrix forms a relatively soft fibrous matrix, it is necessary to handle the material very carefully when performing media changes. **Do not use vacuum aspiration** to remove media from above the hydrogel. Avoid direct contact with the hydrogel.
  - When adding media and/or other components to the sample, place the end of the pipette tip towards the top of the culture vessel wall.
  - When using culture vessels with large surface areas (e.g., ≥35 mm dishes), perform fluid handling steps with extreme care to avoid disruption of the hydrogel.

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences

## CORNING

### **GENERAL USE**

- See table on page 15 to determine appropriate volume for desired Permeable Support (Cell Culture Insert) or Multiple Well Plate format.
- Corning<sup>®</sup> PuraMatrix<sup>™</sup> can be used in the undiluted form at a concentration of 1% (w/v) or diluted as described in **PROCEDURES FOR USE**. Concentrations below 0.5% are suitable for many applications. For example, we recommend 0.25% for surface plating of most cell types, and a final concentration of 0.15% for encapsulating neurons and endothelial cells.
- To supplement the hydrogel with ECMs, growth factors, and/or other bioactive molecules, refer to **Part C (ECM Supplementation Protocols)**.
- If using serum-free media for culturing cells, it may be necessary to deactivate the trypsin following cell isolation by trypsinization.

IMPORTANT: To achieve optimal cell growth and differentiation, it is necessary to determine the appropriate mixture of Corning PuraMatrix Peptide Hydrogel and bioactive molecules (e.g., growth factors, ECM proteins, and/or other molecules). In general, it may be preferable to add growth factors to the media rather than encapsulate them in the hydrogel for most applications.

### PROCEDURES FOR USE

### A) Surface Plating of Adherent Cells in Cell Culture Plates

**Note:** Many cell types that are plated on the surface of Corning PuraMatrix can assume the appropriate and expected 3D structures (e.g., hepatoycyte spheroids, embryonic stem cell colonies and embryoid bodies, neuronal cell will branch). The 2D surface plating method is an ideal tool for screening both morphological and functional characteristics of any given cell type.

### Corning PuraMatrix Preparation and Gel Formation

- Decrease the viscosity of the Corning PuraMatrix stock solution (1% w/v) by vortexing or by sonication in a bath sonicator as described above. Prepare the appropriate volume of Corning PuraMatrix in a microtube by diluting the stock with sterile dH<sub>2</sub>O.
- 2. Prepare a sufficient volume of the desired concentration of Corning PuraMatrix (e.g., 250  $\mu$ l for a 24-well plate or 50  $\mu$ l for a 96-well plate) by dilution with sterile dH<sub>2</sub>O. Add to the surface of the well and promote gelation by carefully and slowly adding medium to each well (e.g., 500  $\mu$ l for a 24-well plate or 100  $\mu$ l for a 96-well plate).
- Place plate in an incubator for 30-60 minutes to complete the gelation of the Corning PuraMatrix. For higher percentage hydrogels, 30 minutes is sufficient; for 0.15-0.25% hydrogels, a 1-hour equilibration is preferable. After the hydrogel has assembled, carefully change the medium as described above. Change the medium two times over a period of 1 hour

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences

CORNING

to

equilibrate the growth environment to physiological pH. If necessary, the equilibrated sample can be stored overnight at 37°C after the final media change.

### Cell Plating

Trypsinize cells and spin down desired number (typically  $4-16 \times 10^4$  cells/cm<sup>2</sup> final concentration for most cell types). Resuspend cells in tissue culture media at an appropriate concentration of cells in a final volume of 200-350 µl. Carefully add the cell suspension to the top of the hydrogel.

### B) Surface Plating of Adherent Cells in Cell Culture Inserts

The following protocol is based on the use of a Falcon<sup>®</sup> 24-well Cell Culture Insert (Cat. No.353095) and a Falcon Cell Culture Insert Companion Plate (Cat. No. 353504).

### **Corning PuraMatrix Preparation and Gel Formation**

- Decrease the viscosity of the Corning<sup>®</sup> PuraMatrix<sup>™</sup> stock solution (1% w/v) by vortexing or by sonication in a bath sonicator as described above. Prepare the appropriate volume of Corning PuraMatrix in a microtube by diluting the stock with sterile dH<sub>2</sub>O. We recommend a total volume of 100 µl per cell culture insert at a final concentration of 0.25% for most cell types. Mix by gentle pipetting.
- Place desired number of cell culture inserts in the insert companion plate. Add 250 µl of media to the lower chambers of each 24-well insert by pipetting down the side of the outer well between the insert and the well. Pipet carefully to avoid bubbles under the insert. The media should just touch the bottom of the insert. See Table on page 15 for cell culture insert volume recommendations.
- 3. Pipet 100 µl of diluted Corning PuraMatrix into the center of each of the inserts. Wait at least 30 minutes for the Corning PuraMatrix to gel.
- VERY CAREFULLY layer 400 μl media onto the gel with a pipet tip (use a P200 pipetter and add 200 μl (two times) for smaller droplet size) by running it slowly down into the inner wall of the insert.
- 5. Follow steps 10-13 in **Part E** to complete pH equilibration of the Corning PuraMatrix prior to cell plating.

### **Cell Plating**

- Trypsinize cells and spin down desired number (typically 4-16x10<sup>4</sup> cells/cm<sup>2</sup> final concentration for most cell types). Resuspend cells in tissue culture media at an appropriate concentration of cells in a final volume of 200-350 μl.
- 2. Remove the media above the hydrogel and from the well as described above.

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences

For Research Use Only. Not for use in diagnostic or therapeutic procedures. For a listing of trademarks, visit www.corning.com/lifesciences/trademarks © 2013 Corning Incorporated

- 3. GENTLY pipet the cells into the inserts using a P200 pipetter, being careful not to disturb the hydrogel.
- 4. Add 700-900 µl of media below each insert.
- 5. For further cell feeding, change the media VERY GENTLY approximately every two days as described above.

### C) ECM Preparation for Supplementation of Corning<sup>®</sup> PuraMatrix<sup>™</sup> Peptide Hydrogel

**Note:** Prior to use, decrease the viscosity of the Corning PuraMatrix stock solution (1% w/v) by vortexing or by sonication in a bath sonicator for 30 minutes every time the stock is used. If air bubbles are present, centrifuge an aliquot of the stock solution in a 1.5 ml microtube.

### ECM Supplementation for Surface Plating of Cells

**Fibronectin Protocol:** (Since fibronectin is not soluble in the absence of salt, it must be "coated" onto the fibers).

- 1. Resuspend 5 mg Corning Fibronectin (Cat. No. 356008) in 5 ml sterile  $dH_2O$ .
- Dialyze overnight at Room Temperature (RT) using 10,000 MW cutoff Slide-A-Lyzer® (Pierce Cat. No. 66453), against dPBS without Ca<sup>2+</sup> and Mg<sup>2+</sup>.
- 3. Prior to addition of fibronectin, mix Corning PuraMatrix thoroughly with gentle pipetting and add to well of appropriate plate at desired concentration. See Table on page 15 for suggested volumes.
- 4. Gel Corning PuraMatrix by adding Fibronectin/PBS solution (1 mg/ml) on top of the Corning PuraMatrix for 60 minutes.
- 5. Carefully rinse with cell culture media for 5 minutes, aspirate and then carefully add cells to the top of the matrix.

**Collagen I Protocol:** (Collagen I is soluble in 1 mM HCl and thus can be incorporated into the Corning PuraMatrix before gelation):

- 1. Dialyze collagen I (stock solution in acetic acid, Cat. No. 354236) using 10,000 MW cutoff Slide-A-Lyzer® against 1 mM HCl overnight at 4°C.
- 2. Add the required amount of collagen I/HCI to Corning PuraMatrix and mix with gentle pipetting.
- 3. Add to well.
- 4. Gel with cell culture media by adding media carefully to side wall of cell culture insert or well of plate.

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences

## CORNING

5. Once gel has formed, remove media carefully with a pipet and change two more times over the next one hour to equilibrate the Corning PuraMatrix to physiological pH prior to plating the cells on top of the matrix.

Laminin Protocol: (Laminin is provided in 0.05 M Tris-HCl, 0.15 M NaCl, pH 7.4).

- 1. Low concentration laminin [Laminin, mouse (Cat. No. 354232) and Ultra-pure Laminin, mouse (Cat. No. 354239)].
  - a. Mix with cell culture media at the desired concentration.

b. Carefully add the mixture directly on top of Corning<sup>®</sup> PuraMatrix<sup>™</sup> (prepared in the presence or absence of cells) in the cell culture insert or plate to promote gelation.

2. Laminin/Entactin, High Concentration (HC), 10 mg/ml (Cat. No. 354259).

Add required volume of HC Laminin directly to Corning PuraMatrix. At dilutions of HC Laminin that result in a final concentration of  $\leq$ 1 mM salt, gelation of Corning PuraMatrix will not occur. To induce gelation after the addition of Laminin, add cell culture media and equilibrate the pH as described in the collagen I protocol.

### ECM Supplementation of Corning PuraMatrix for Cell Encapsulation

# Detailed instructions for cell encapsulation in Corning PuraMatrix can be found in Parts D and E. These should be followed, with the modifications detailed below, to introduce ECMs into the cell encapsulation mixture.

For Corning PuraMatrix concentrations ≤0.25%, a cell culture insert is recommended to facilitate media changes and ensure that the hydrogel is not disrupted. Corning PuraMatrix concentrations below 0.25% are recommended for the encapsulation of neurons and endothelial cells. With other cell types, testing a range of hydrogel concentrations for encapsulation is recommended.

### To add ECMs to cell encapsulation mixture for cell culture plates or cell culture inserts:

- 1. Prepare ECMs as described above.
- Add laminin or collagen (fibronectin is not suitable for encapsulation, but it can be coated onto samples of cells that are encapsulated within Corning PuraMatrix) to the cell/sucrose mixture to a final concentration of 1-5 μg/ml (2X desired final concentration).
- 3. For cell encapsulation, follow protocols in **Part D** (cell culture plates) or **Part E** (cell culture inserts).

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences



## D) 3D Cell Encapsulation in Cell Culture Plates (with or without ECM protein supplementation)

- 1. Decrease the viscosity of the Corning<sup>®</sup> PuraMatrix<sup>™</sup> stock solution (1% w/v) by vortexing or by sonication in a bath sonicator for 30 minutes every time the stock is used. If air bubbles are present, centrifuge an aliquot of the stock solution in a 1.5 ml microtube.
- 2. Prepare the appropriate volume of desired concentration of Corning PuraMatrix in a microtube by diluting the stock with sterile 20% sucrose to generate a 2X concentration of Corning PuraMatrix in 10% sucrose. See Table on page 15 for cell culture plate volume recommendations.
- 3. Trypsinize cells and spin down desired number (typically 5x10<sup>5</sup>-1x10<sup>6</sup> cells/ml final concentration for most cell types).
- 4. Remove media from cell pellet and resuspend cells to desired density in sterile 10% sucrose.
- 5. Collect cells again by centrifugation and resuspend the cells in 10% sucrose at 2X the final desired cell concentration.
- 6. Add equal volumes of 2X Corning PuraMatrix and 2X cell/sucrose mixture, mix, and then add to the center of the well carefully, without introducing bubbles.
- 7. Initiate gelation of the Corning PuraMatrix by gently running culture media down the side of the well on top of the hydrogel.
- 8. Repeat until cells have been plated in all wells.
- 9. VERY GENTLY change the media two times over the next one hour to further equilibrate the pH of the hydrogel. DO NOT USE a vacuum aspirator and only remove about 2/3 to 3/4 of the media to avoid disrupting the Corning PuraMatrix.
- 10. Cells can be fed as necessary using the procedure described in step 9.

## E) 3D Cell Encapsulation in Cell Culture Inserts (with or without ECM protein supplementation)

The following protocol is based on the use of a Falcon<sup>®</sup> 24-well Cell Culture Insert (Cat. No. 353095) and a Falcon Cell Culture Insert Companion Plate (Cat. No. 353504).

For Corning PuraMatrix concentrations ≤0.25%, a cell culture insert is recommended to facilitate media changes and ensure that the hydrogel is not disrupted. Corning PuraMatrix concentrations below 0.25% are recommended for the encapsulation of neurons and endothelial cells. With other cell types, testing a range of hydrogel concentrations for encapsulation is recommended.

### CORNING

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences

### **Corning PuraMatrix Preparation**

- 1. Decrease the viscosity of the Corning PuraMatrix stock solution (1% w/v) by vortexing or by sonication in a bath sonicator for 30 minutes every time the stock is used. If air bubbles are present, centrifuge an aliquot of the stock solution in a 1.5 ml microtube.
- Prepare the appropriate volume of the desired concentration of Corning<sup>®</sup> PuraMatrix<sup>™</sup> in a microtube by diluting the stock with sterile 20% sucrose to generate a 2X concentration of Corning PuraMatrix in 10% sucrose. We recommend a total volume of 100 µl per 24-well cell culture insert (i.e., 50 µl of 2X Corning PuraMatrix in 10% sucrose and 50 µl of cells at 2X in 10% sucrose). See Table on page 16 for cell culture insert volume recommendations.
  - a. For example, to encapsulate neurons and endothelial cells in a final concentration of 0.15% Corning PuraMatrix, prepare a stock of 0.30% Corning PuraMatrix in 10% sucrose (for 1 ml, combine 500 µl of 20% sucrose, 300 µl Corning PuraMatrix, and 200 µl of sterile H<sub>2</sub>O).
- 3. Mix by gentle pipetting.
- 4. We recommend setting up a separate tube to mix the cells with Corning PuraMatrix for each cell culture insert or well. The desired volume of Corning PuraMatrix can be distributed to all the tubes at once.

**Cell Seeding** – Since Corning PuraMatrix exhibits low pH, work quickly to minimize the amount of time that cells are in contact with the material prior to the addition of media. When plating multiple wells, we recommend that each microtube is prepared individually, as opposed to in a group, to avoid long exposure of the cells to the Corning PuraMatrix in the absence of media.

- 1. Trypsinize cells and spin down desired number (typically 5x10<sup>5</sup>-1x10<sup>6</sup> cells/ml final concentration for most cell types).
- 2. Place desired number of cell culture inserts in the insert companion plate.
- 3. Add 250 µl of media to the lower chambers of each insert by pipetting down the side of the outer well between the insert and the well. Pipet carefully to avoid bubbles under the insert. The media should just touch the bottom of the insert (**Figure 1**).
- 4. Remove media from cell pellet and resuspend cells to desired density in sterile 10% sucrose.
- 5. Collect cells again by centrifugation and resuspend the cells in 10% sucrose at 2X the final desired cell concentration.
- 6. QUICKLY mix diluted 2X Corning PuraMatrix in 10% sucrose in a microtube with an equal volume of 2X cell suspension in 10% sucrose using gentle pipetting. Do not vortex or invert once cells are added.

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences

For Research Use Only. Not for use in diagnostic or therapeutic procedures. For a listing of trademarks, visit www.corning.com/lifesciences/trademarks © 2013 Corning Incorporated

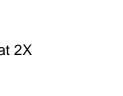



Figure 1




7. QUICKLY pipet Corning PuraMatrix/cell mixture into a 24-well insert (**Figure 2**). When dispensed into the center of the insert, the material will spread evenly over the surface.

Figure 2

- 8. Wait five minutes for the Corning<sup>®</sup> PuraMatrix<sup>™</sup> to settle in the well. During this period, set up the next insert as described in steps 6-7.
- Return to the first insert and VERY CAREFULLY layer 400 μl media onto the gel with a pipet tip (use a P200 pipetter and add 200 μl (two times) for smaller droplet size) by running it slowly down the inner wall of the insert (Figure 3).

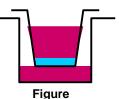

### Flowchart of steps 6 through 9

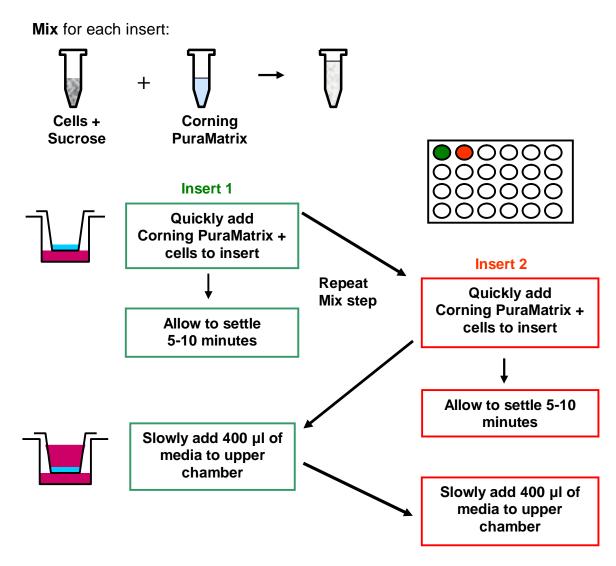

Repeat until all inserts are seeded, then continue with step 10 of the protocol.

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences

CORNING

For a listing of trademarks, visit www.corning.com/lifesciences/trademarks


© 2013 Corning Incorporated

For Research Use Only. Not for use in diagnostic or therapeutic procedures.

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences



10. When all inserts have been filled, place the plate in an incubator for 30-60 minutes to complete the gelation of the Corning<sup>®</sup> PuraMatrix<sup>™</sup>. For higher percentage hydrogels, 30 minutes is sufficient; for 0.15-0.25% hydrogels, a one-hour equilibration is preferable. The media levels may equilibrate to some extent between the insert and the well during this period (**Figure 4**).

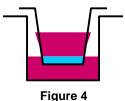

Figure

11. Change the media without disturbing the hydrogel. As the gelled Corning PuraMatrix is difficult to see, this is best accomplished by removing approximately 2/3 of the media above the hydrogel with a P200 pipetter and keeping the tip above the bottom of the insert (Figure 5). DO NOT USE A VACUUM to aspirate media from above the hydrogel as it will disturb the Corning PuraMatrix scaffold and cells.

The media in the wells below the insert can be removed with a vacuum aspirator or a pipet tip.

- 12. VERY GENTLY replace the media removed from the inserts, and also in the wells, as described above. Place the plate in the incubator for 30 minutes and then repeat the media change as described in step 11 and incubate again. This will equilibrate the pH of the Corning PuraMatrix and remove excess sucrose. After the second change, place 200-350 µl of media into the insert and 700-900 µl into the well.
- 13. For further cell feeding, change the media VERY GENTLY approximately every two days as described above.

### F) Cell Drop Culture Protocol

This technique, which allows the formation of small beads containing cells encapsulated within Corning PuraMatrix, can be used as a quick screen for optimizing the conditions for encapsulating the cell type of interest within Corning PuraMatrix (in the presence or absence of ECM proteins). The cell drop culture protocol can also be used to quickly assess cell viability, cell morphology, and/or immunostaining of cells encapsulated within Corning PuraMatrix.

- 1. Decrease the viscosity of the Corning PuraMatrix stock solution (1% w/v) by vortexing or by sonication in a bath sonicator for 30 minutes every time the stock is used. If air bubbles are present, centrifuge an aliquot of the stock solution in a 1.5 ml microtube.
- 2. Prepare the appropriate volume of Corning PuraMatrix in a microtube by diluting the stock with sterile 20% sucrose to generate a solution with 2X the desired final concentration of Corning PuraMatrix in 10% sucrose. Mix by gentle pipetting.
- 3. Grow cells to 75-80% confluence and harvest using standard procedures.

## CORNING

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences

- If ECM proteins are to be added to the drops, prepare as described in Part C (ECM Supplementation Protocol). The final salt concentration should be ≤1 mM to prevent gelation of the Corning<sup>®</sup> PuraMatrix<sup>™</sup>.
- Gently pellet the cells at 1000 x g, wash with 5 ml of 10% sucrose (filter sterilized, 0.2 μm membrane), re-pellet, and then resuspend with 10% sucrose to 2X desired final concentration (generally 5x10<sup>5</sup> 1x10<sup>6</sup> cells/ml). Mix ECM proteins to 2X final concentration at this step. Use gentle pipetting to mix the sample.
- 6. Combine the cell/sucrose suspension with an equal volume of Corning PuraMatrix and mix with gentle pipetting. Dispense 15 μl of the cell/hydrogel mixture along the side of the well on the surface of the media (1 ml of complete medium), which initiates the assembly of the hydrogel. Multiple "drops" can be dispensed in the well containing medium.
- 7. To assess cell viability after the desired time in culture, the drop culture can be analyzed using the LIVE/DEAD Viability/Cytotoxicity Assay Kit (Molecular Probes Cat. No. L-3224). Collect the cell/gel particles with a spatula from the surface of the media. Place the bead sample on a glass slide, and then add 50 µl of LIVE/DEAD Calcein AM and EthD-1 solution.
- 8. Incubate the slide for 30 minutes at room temperature and then examine using an imaging microscope for the presence of green (LIVE) and red (DEAD) cells.

### G) Cell Recovery for Sub-Culturing or Biochemical Analyses

- 1. Mechanically disrupt Corning PuraMatrix and cells in the well or cell culture insert by repeated pipetting of the media and gel.
- 2. Transfer to a 15 ml conical tube. Use another half to full volume of PBS to rinse out the well or cell culture insert and combine in the tube.
- 3. Centrifuge at low speed for five minutes. Discard supernatant. The pellet at the bottom of the tube contains cells and Corning PuraMatrix fragments.
- 4. Resuspend pellet in 2 ml PBS, spin, and collect pellet again.
- 5. Resuspend pellet in 1 ml of trypsin-EDTA and incubate at 37°C for 5-10 minutes. This will help separate cells that are still attached to each other.
- 6. Add 5 ml media and repellet.
- Add 500 µl media (for replating) or PBS (for biochemical analyses) to the pellet and resuspend by pipetting. Centrifuge and discard supernatant. At this point, some cells may have residual "tags" of Corning PuraMatrix attached which cannot be removed.

## CORNING

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences

- 8. For biochemical and molecular analyses, add chosen buffer to pellet for cell solubilization.
- For cell plating, transfer the pellet in media to a new cell culture plate or insert. Cells can be reencapsulated (see **Parts D and E**), surface plated on a layer of Corning<sup>®</sup> PuraMatrix<sup>™</sup> (see **Parts A and B**) or seeded directly onto a plastic plate (uncoated or coated with an ECM protein).

### H) Fluorescent Staining of Samples Cultured with Corning PuraMatrix Peptide Hydrogel

### **IMMUNOSTAINING**

This protocol can be used with 3D or 2D surface cultures of cells in inserts or plates. For certain cell culture formats (e.g., 24- and 96-well microplates with a clear bottom), it may be possible to directly image the samples without transferring the gel to a glass slide. In some cases, it may be necessary to transfer the gel to a glass side for imaging analysis [e.g., when using a cell culture insert with translucent membrane (e.g., Cat. No. 353495)].

**Note:** To achieve a strong and specific signal, it is necessary to thoroughly block non-specific binding of antibodies to Corning PuraMatrix as well as allowing diffusion of antibodies through the peptide scaffold. Therefore, extended blocking and incubation times, higher antibody concentrations, and multiple washes are required.

- 1. Gently remove media from culture well.
- 2. Fix cells in gels in 4% paraformaldehyde for 30 minutes.
- 3. Wash in PBS.
- 4. Incubate in Block Solution (PBS + 10% Fetal Bovine Serum) for 12-16 hours.
  - a. Block solution should be changed three to four times every few hours and then left overnight. Additional changes can be done if necessary.
- 5. Add primary antibody in block solution. Incubate overnight at 4°C.
  - a. A higher concentration of primary antibody is recommended for cells present in a 3D Corning PuraMatrix sample compared to that used for cells in 2D culture. Determine the optimal concentration by titrating the antibody.
  - b. Process one sample **without primary antibody** as a negative control to assess background staining.
  - c. Use sufficient volume to completely soak the 3D gel.

### CORNING

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences

- 6. Wash in block solution at least four times (two hours per wash).
- 7. Add secondary antibody in block solution.
  - a. As described in step 5a above, a higher concentration of secondary antibody is recommended for cells in Corning<sup>®</sup> PuraMatrix<sup>™</sup>.
  - b. Incubate for four hours.
- 8. Wash again in PBS or block solution multiple times (four to six times) for at least one hour per wash.
- 9. Perform imaging analysis of the samples in the cell culture vessel directly. If necessary, carefully **remove the gel** from the vessel with spatula and place on a glass slide for microscopic analysis. If using a low percentage gel (e.g., <0.5%), use extreme care when transferring the gel.

### MORPHOLOGICAL STAINING: RHODAMINE-PHALLOIDIN/DAPI

- 1. Gently remove media from insert or culture well.
- 2. Fix cells in gels in 4% paraformaldehyde for 20 minutes.
- 3. Gently remove fixative and discard.
- 4. Wash 1X with PBS
- 5. Add desired stain.
  - a. We recommend 160 nM rhodamine-phalloidin in PBS, 300 µl per well or per insert.
  - b. Cover and leave overnight.
  - c. Add 1X DAPI for one hour and then rinse 2X with PBS.
- 6. Remove PBS and perform imaging analysis as described above for immunostaining.

### I) Tips for *In Vivo* Delivery of Cells and/or Bioactive Molecules

- Corning PuraMatrix can be readily handled using small or large bore needles and catheters. Small-bore needles may generate bubbles. To avoid the introduction of bubbles *in vivo*, use extreme care when filling and injecting samples with needles that are smaller than 20G.
- The material can be re-filtered with a 0.20 or 0.45 µm membrane filter prior to use.
- The material should be vortexed or sonicated prior to use to reduce viscosity.

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences

## CORNING

While Corning<sup>®</sup> PuraMatrix<sup>™</sup> can be injected as provided at low pH, some injection sites may not tolerate exposure to an acidic sample. Also, cells are sensitive to exposure to low pH. Therefore, the following protocol can be used to buffer and gel the material prior to injection into animals.

To prepare a 200  $\mu I$  sample for injection, begin by aliquoting 50  $\mu I$  of a cell/10% sucrose suspension into a small tube.

- 1. Add 50 µl of Corning PuraMatrix and mix gently.
- Immediately layer 100 μl of PBS (containing Ca<sup>2+</sup> and Mg<sup>2+</sup>) carefully over top of cell/hydrogel mixture.
- 3. Allow to sit at room temperature for 1 minute. Mix gently with a P200 pipet tip.
- 4. Using a 20-22 g needle, pull entire mixture into 1 cc syringe. Avoid air bubbles by slowly and carefully collecting the liquid.
- 5. A needle change is recommended for injection. In general, use a 30 g or 26 g needle for tissue or subcutaneous injections, respectively.

### Notes:

- The volumes can be scaled according to requirements of the animal used and the injection site.
- The number of cells injected should be determined empirically according to the goals of the experiment.

### **Recommended Plating Volumes for Corning PuraMatrix Peptide Hydrogel**

	Growth Area (cm <sup>2</sup> )	Volume Per Well (µl)				
Multiple Well Plates						
6-well	9.60	1200				
24-well	2.00	250				
96-well	0.32	50				
	Growth Area (cm <sup>2</sup> )	Volume Per Insert (µl)				
Permeable Supports (Cell Culture Inserts)						
6-well	4.20	700				
24-well	0.30	50				
96-well	0.08	14				

Note: The final volumes of Corning PuraMatrix for plating (Corning PuraMatrix/cell/sucrose mixture) should be doubled for cell encapsulation.

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences



For Research Use Only. Not for use in diagnostic or therapeutic procedures.

For a listing of trademarks, visit www.corning.com/lifesciences/trademarks © 2013 Corning Incorporated

### STORAGE

Stable for when stored at 4-30°C.

To place an order in the U.S., contact Customer Service at: tel: 800.492.1110, fax: 978.442.2476; email: <u>CLSCustServ@corning.com</u>.

For technical assistance, contact Technical Support at: tel: 800.492.1110, fax: 978.442.2476; email: <u>CLSTechServ@corning.com</u>.

Outside the U.S., contact your local distributor or visit <u>www.corning.com/lifesciences</u> to locate your nearest Corning office.

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences



### **Revision History**

	ECO Number			
Rev		Description	Revised By	Date
1.0	ECO12508	Initial Issue	Paula Flaherty	05/03/04
2.0	ECO18731	Added in in vivo injection info, clarified wording	Paula Flaherty	08/26/04
3.0	ECO47846	Technical information, graphics, and format	Marshall	5/31/06
		changes	Kosovsky/Lee	
			Stowell	
4.0	ECO81910	Technical revisions, process clarification under	S. Bunch/M.	4/9/08
		procedures Part D, added reference on pg. 2	Kosovsky	
5.0	ECO119495	Updated to new logo, updated tech support and	Jen Amaral	5/20/10
		customer service info		

Discovery Labware, Inc., Two Oak Park, Bedford, MA 01730, Tel: 1.978.442.2200 (U.S.) CLSTechServ@Corning.com www.corning.com/lifesciences

