


Lonza GelBond® Film

Lonza

A transparent, flexible polyester film designed to support agarose gels. Gels cast on GelBond® film remain permanently attached to the film through electrophoresis or immunodiffusion and all subsequent fixing, staining, destaining, and drying procedures. GelBond® film is available either as pre-cut sheets or rolls. Note: polyester films will not transmit light of less than 310nm, and will fluoresce at higher wavelengths.

- Agarose gels cast on GelBond® film retain their original dimensions during staining and after drying
- Gels, particularly thin ones, are easier to handle during staining, destaining, and drying when supported
- Gel orientation can be recorded directly on the GelBond® film prior to casting

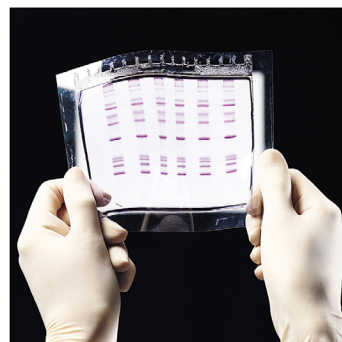


Code	Dims, w x l, mm	Pack	Price
LZ53734	85 x 100	100 sheets	£376.00
LZ53745	110 x 125	100 sheets	£548.00
LZ53746	100 x 150	100 sheets	£615.00
LZ53748	110 x 205	100 sheets	£829.00
LZ53749	160 x 180	100 sheets	£857.00
LZ53759	125 x 245	100 sheets	£926.00
LZ53761	124 x 258	100 sheets	£978.00
LZ53740	102 x 16500	Roll	£559.00
LZ53750	127 x 16500	Roll	£720.00
LZ53780	203 x 16500	Roll	£957.00

Lonza GelBond® PAG Film

Lonza

A transparent, flexible polyester film designed to support polyacrylamide or MDE™ gels. The acrylamide monomers covalently attach to the coating on the film during the polymerisation reaction. Gels remain permanently attached to the film through electrophoresis and all subsequent fixing, staining, destaining, and drying procedures.

- Polyacrylamide gels retain their original dimensions during staining and after drying
- Gels, particularly thin ones, are easier to handle during staining, destaining, and drying when supported
- Gel orientation can be recorded directly on the GelBond® PAG film prior to casting

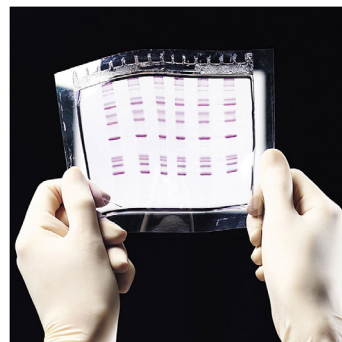

Code	Dims, w x l, mm	Pack	Price
LZ54711	138 x 158	50 sheets	£309.00
LZ54723	160 x 180	50 sheets	£355.00
LZ54727	124 x 258	50 sheets	£432.00
LZ54729	165 x 220	50 sheets	£566.00
LZ54731	199 x 264	50 sheets	£725.00
LZ54733	203 x 260	50 sheets	£610.00
LZ54746	350 x 430	10 sheets	£419.00