


Astell Compact Autoclaves



Front loading autoclave is ideal for traditional laboratory sterilisation procedures, with the option of a manual or autofill version. Programmes can be set to accommodate all common laboratory sterilisation requirements, with time and temperature easily adjustable.

All models feature:

- Multi-programmable touch screen controller (up to 50 cycles)
- Temperature range 100 to 138°C (0.2 to 2.4 bar)
- 316 stainless steel pressure vessels
- Fast action door mechanism
- Designed for maximum chamber capacity
- No plumbing, and operate from a standard single phase power supply
- Front loading, cylindrical chambers

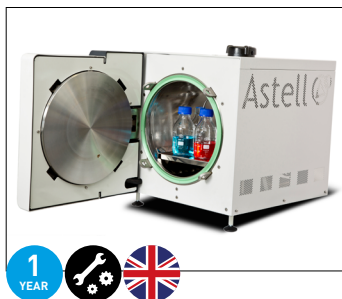

Classic Version

This 'Classic' model requires manual water-fill with condensate leaving the vessel via a silicone tube that can be connected to an optional heat-resistant bottle. Media 'Holdwarm' is a standard feature on this version.

Code	Alt Ref	Capacity, L	Chamber Dims, Ø x d, mm	Overall Dims, w x d x h, mm	Weight, kg	Price
AUT4400	AMB420	33	346 x 355	520 x 805 x 600	96	£8830.00
AUT4402	AMB430	43	346 x 465	520 x 805 x 600	96	£10235.00
AUT4404	AMB440	63	346 x 668	520 x 1005 x 600	105	£12065.00



Ecofill Version

The Ecofill version of Astell's Benchtop range has an integral reservoir which, once filled provides water for up to 15 cycles and acts as a collection vessel for condensate which is automatically returned from the chamber at the end of each cycle. This means that the unit is totally self-contained - all you need is a standard power supply.

Code	Alt Ref	Capacity, L	Chamber Dims, Ø x d, mm	Dims, w x d x h, mm	Weight, kg	Price
AUT4300	AMB220	33	346 x 355	520 x 805 x 600	96	£10016.00
AUT4220	AMB230	43	346 x 355	520 x 805 x 600	96	£11523.00
AUT4304	AMB240	63	346 x 355	520 x 1005 x 600	105	£13380.00