



## Azenta Life Sciences FrameStar® 384 Well Skirted PCR Plates



- FrameStar® 384 well skirted PCR plate PP wells, PC frame, cut corner A2
- Working volume: 30µL
- Total well capacity: 55µL
- Compatible with the majority of 384 well block PCR, qPCR and sequencing instruments
- Available with frosted wells for increased qPCR signal intensities
- Available in different frame colours

Code	Alt Ref	Well Colour	Frame/Arm Colour	Pack	Price
PCR0430	4ti-0384/C	Clear	Clear	50	£218.00




## Azenta Life Sciences Semi-Automated Sheet Heat Sealer



A semi-automatic heat sealer which is compatible with a wide range of seals and plates of differing designs and heights. With variable temperature and time settings, sealing conditions are easily optimised to produce a 100% seal, eliminating sample loss. Plate and seal are placed on the holder, the "Operate" button pressed and the drawer automatically closes. The sealing process is controlled by an electric mechanism; this unique feature eliminates operator variation, giving a consistent seal every time and reducing the damage caused to microplates by over-sealing.

- Rapid heating element to allow for a quick start up within minutes
- Real time temperature display to allow monitoring of heating process
- Automatically switches to standby mode when idle; heating element temperature reduces to 60°C to save energy
- When idle for longer periods, it automatically switches off for added safety
- Triggering of standby and switch off modes can be adjusted to suit user's needs
- Unique and cost saving adapter system
- Compatible with all SBS microplates
- RS-232 port enables unit integration into most automated systems

Code	Temp Range, °C	Dims, w x l x h, mm	Weight, kg	Price
59-2000	60 to 200	181 x 275 x 310	9	£3819.00

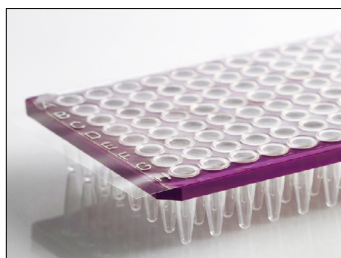

## Azenta Life Sciences Adhesive Sealing



- Durable transparent polyester film with a strong adhesive layer, peelable
- Enables high integrity and efficiently prevents sample evaporation
- Recommended for PCR and other optical applications as well as sample storage qPCR seal
- Clear, non-sticky film with a strong pressure-activated adhesive, peelable
- Recommended for qPCR and other imaging techniques including crystallisation
- Also available as half plate seals to match the size of our FrameStar® 192 well semi-skirted PCR plate

Code	Alt Ref	Description	Dimensions, w x l	Pack	Price
PCR0516	4ti-0500	Adhesive PCR plate seals	135 x 80mm	100 sheets	£114.00
PCR0518	4ti-0500/8	Adhesive PCR plate seals with 12 vertical perforated strips of 8 wells	115 x 100mm	100 sheets	£115.00
PCR0526	4ti-0500/12	Adhesive PCR plate seals with 8 horizontal perforated strips of 12 wells	137 x 71mm	100 sheets	£115.00
PCR0564	4ti-0560	qPCR adhesive clear plate seals sheets	140 x 77mm	100 sheets	£167.00
PCR0566	4ti-0561	qPCR pressure sensitive adhesive seal roll	100m x 80mm	Roll (approx. 700 seals)	£842.00