

# Consumables for the Shandon Cytospin® Cytocentrifuges

## Single CytoSep™ Cytology Funnel



M964-10FW



M964-10FT



Fully compatible with the Shandon Cytospin® Centrifuge, Simport® CytoSep™ Cytology Funnels can be used to deposit a thin layer of cells in a clearly defined area of a microscope slide. The filter card absorbs any excess fluid. These Cytology Funnels have the filter cards pre-attached for consistent, reliable results. No alignment necessary! All disposable Funnels are packaged with closure caps to seal in specimen for added protection. All components also available separately.

For sample volumes of up to 0.5 ml, use the Simport® Single CytoSep™ Cytology Funnel with a White Filter Card and Cap. It provides a cell deposition area of 6 mm (28 mm<sup>2</sup>). Can be use with all stainless steel slide clips.

For samples volumes of up to 0.4 ml, such as Spinal Fluids for example, use the Single CytoSep™ Cytology Funnel with Brown Filter Card and Cap. It allows for a slower absorption of fluids. All the individual components are also available separately.



A more convenient packaging design has been developed for individually wrapped and bulk packaged Cytology Funnels, thereby facilitating their handling.

For IVD use

Cat. #	Description	Qty/Pk	Qty/Cs
M964-10FW	Single Funnel with White Filter & Cap	50	500
M964-10FW1	<b>Individually wrapped</b> Single Funnel with White Filter & Cap	100	500
M964-10FT	Single Funnel with Tan Filter & Cap	50	500
M964-1	Single Funnel only	–	500
M965C	Cap only	50	500
M965FW	White Filter Paper for Single Funnel	200	–
M965FT	Tan Filter Paper for Single Funnel	200	–

## Double CytoSep™ Cytology Funnel



M964-20FW



The Simport® Double CytoSep™ Cytology Funnel with disposable sample chamber allows for two samples to be run simultaneously on a single slide and is ideal for immunohistochemistry work. The Double CytoSep™ Cytology Funnel provides a cell deposition area of 6 mm (28 mm<sup>2</sup>) for sample volumes of up to 0.5 ml. The filter card comes pre-attached. Cap is included. Can be use with all stainless steel slide clips. Each component is also available separately.

For IVD use


Cat. #	Description	Qty/Pk	Qty/Cs
M964-20FW	Double Funnel with White Filter & Cap	50	500
M964-20FW1	<b>Individually wrapped</b> Double Funnel with White Filter & Cap	100	500
M964-1D	Double Funnel only	–	500
M965C	Cap only	50	500
M965FWD	White Filter Paper for Double Funnel	200	–

## M964B CytoSep™ Metal Cytology Funnel Clip

Made of Stainless Steel

The Simport® CytoSep™ Metal Cytology Funnel Clip is similar in design to the model previously provided by Shandon. It will hold the Cytology Funnel disposable sample chamber against the microscope slide in the Thermo Scientific Cytospin™ Cytocentrifuge. This Simport® reusable metal clip is autoclavable and will not rust or corrode during decontamination.


Used with all M964  
CytoSep™ Cytology  
Funnels

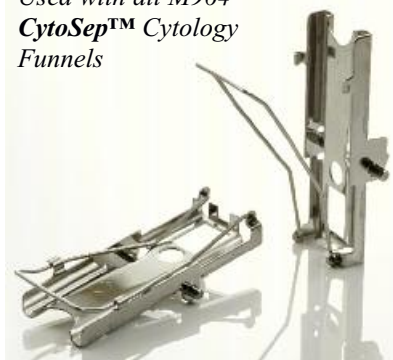

To open, push clip sideways to the right



Insert microscope slide

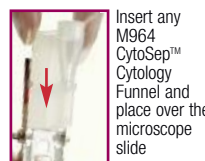

Insert any M964 CytoSep™ Cytology Funnel and place over the microscope slide



Secure clip by pushing it sideways to the left into the hooks

For IVD use

Cat. #	Description	Qty/Pk	Qty/Cs
M964B	Metal Cytology Funnel Clip	–	6