


The dual purpose scoops by Heathrow Scientific are a unique economical solution that combines the function of a scoop, weigh boat, and funnel into one product.

The weighing process has been streamlined, now just sample or retrieve the material in the scoop and place it directly on the scale. This eliminates the need for a weigh boat or weighing paper and saves valuable space in the lab.

- Made from anti-static polystyrene that is waterproof and reduces the chance of materials sticking to the scoop
- Ships and stores flat and when ready to use, is easy to assemble by inserting the locking tab into the hole
- Locking tab feature enables a scoop shape to help retain materials
- Coated to stand up to powders, gels, and semi-solids
- Broad tip design keeps the scoop stable while on the scale and enables materials to be easily scooped and dispensed
- Certificate of Sterility available upon request for sterile scoops
 - Electron beam sterilized products are ideally suited for biologically sensitive applications
 - The dosage requirements for electron beam sterilization that results in sterilized products is determined by procedures that are compliant with ISO 11137-3:2017
- Disposable



Item No.	Color	Size	L x W in	L x W cm	Sterile	UOM	
120656	White	<input type="checkbox"/>	Large	5 x 3.5	12.7 x 9	✓	100 ea/pk, 50 pks/cs
120679	Black	<input checked="" type="checkbox"/>					10 ea/bag, 10 bags/pk, 25 pks/cs
120657	White	<input type="checkbox"/>	Small	5 x 2.5	12.7 x 6.3	✓	100 ea/pk, 50 pks/cs
120680	Black	<input checked="" type="checkbox"/>					10 ea/bag, 10 bags/pk, 25 pks/cs

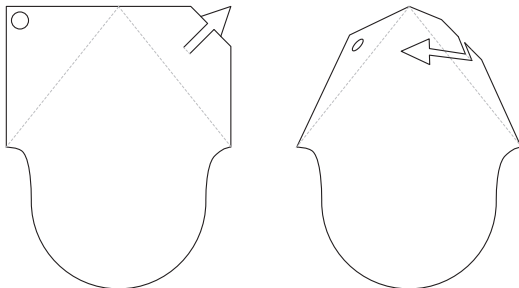

Easily assemble by pulling the arrow tab through the hole.



Sterile Scoops

- Electron beam sterilized
- Lot number traceability for quality assurance