CORNING	Life Sciences	Document Number: LSR00061	Rev.: 3
Description (Class): Reference, Customer Technical Data Sheet Cell Culture Flask, 25 cm ² canted neck, plug seal cap. Total volume – 50 mL			

APPLICATIONS:

- Cell culture
- Sample or reagent storage

FEATURES:

Polystyrene construction

- Optically clear for easy viewing of sample or cells.
- Stable from 4°C to 60°C
- Available in Standard Tissue Culture (TC) (Cat. No. 353014) or Corning[®] Primaria[™] (Cat. No. 353813)
- Frosted writing patch for easy labeling

Two-position plug seal cap

- Biologically inert, high density polyethylene provides a chemically resistant surface
- Cap can be loosely closed to allow gas exchange
- Further tightening snaps cap into sealed position that is liquid-tight
- A reference bar on the outside of the cap can be aligned with flask markings for precise incubation conditions

Molded graduations along both sides of flasks

- Graduated at 5 mL increments to 40 mL
- Allows quick volume estimation

Prominent stacking rims along perimeter of upper surface

- Permits easy interlocked stacking in incubator or culture hood
- Flasks will not slip off one another while being transported

Tissue culture-friendly packaging

- Convenient reseal tab on the bag
- Innovative bag materials that will not scuff or scratch the flask surface
- Double-wall bags provide increased sterility assurance
- Medical-style peel-open bags guarantee that sterility is maintained as flasks are removed from the bag
- Recyclable bag material (low density polyethylene)

STERILITY:

Product is gamma irradiated.

REGULATORY COMPLIANCE:

Falcon[®] products are manufactured under the current ISO 9001 and ISO 13485 Standards, and the current FDA Quality System Regulation 21 CFR Section 820.

RESTRICTIONS:

 Refer to "Thermoplastics Properties Chart," available at: http://catalog2.corning.com/Lifesciences/media/pdf/an_DL_226_Falcon_ThermoPlastics_Broch.pdf


CORNING	Life Sciences	Document Number: LSR00061	Rev.: 3
Description (Class): Reference, Customer Technical Data Sheet Cell Culture Flask, 25 cm ² canted neck, plug seal cap. Total volume – 50 mL			


Material: Polystyrene

Related Catalog Numbers:

Cat. No.	Packout Description	Сар	Neck Style	Graduations	Surface
353014	20/bag, 200/case	Plug Seal	Canted	Numeric/molded	Stand TC
353813	20/bag, 200/case	Plug Seal	Canted	Numeric/molded	Corning [®] Primaria [™] TC

Dimensions = mm (inches)

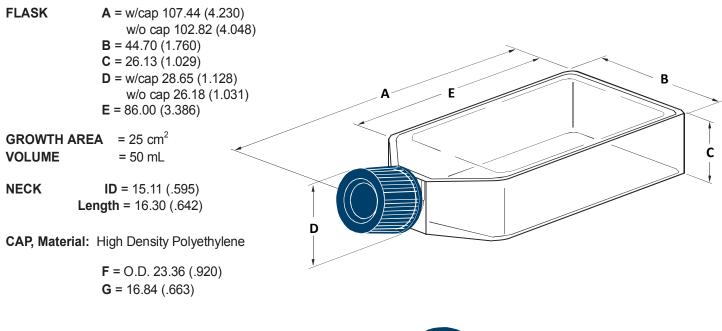

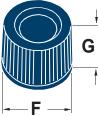

Warranty/Disclaimer: Unless otherwise specified, all products are for research use only. Not intended for use in diagnostic or therapeutic procedures.

For a listing of trademarks, visit www.corning.com/lifesciences/trademarks. © 2013 Corning Incorporated