

# Model BF 720 | Standard-Incubators with forced convection

The BINDER incubator of the BF Avantgarde. Line series with forced convection is suitable for all gentle incubation applications, particularly under a full load and large batch throughputs. This incubator is extremely homogeneous and has quick recovery times.

## **BENEFITS**

- · Excellent temporal and spatial temperature accuracy
- · High process reliability
- · 100 °C disinfection routine



Model 720

## MAIN FEATURES

- Temperature range: ambient temperature plus 10  $^{\circ}\text{C}$  to 100  $^{\circ}\text{C}$
- High temperature accuracy thanks to APT.line<sup>™</sup> technology
- · Forced convection
- · Controller with LCD display
- · Electromechanical control of the exhaust air flap
- Inner door made of tempered safety glass

- · 2 chrome-plated racks
- 4 stable castors, two with brakes, from 720 liters
- Class 3.1 integrated independent temperature safety device (DIN 12880) with visual alarm
- Ergonomic handle design
- USB port for recording data

#### ORDERING INFORMATION

Interior volume [L]	Voltage	Option model	Version	ArtNo.
734	230 V 1~ ph 50/60 Hz	Standard	BF720-230V	9010-0321
1 34	240 V 1~ ph 60 Hz	Standard	BF720UL-240V	9010-0322

## TECHNICAL DATA

Description	BF720-230V <sup>1</sup>	BF720UL-240V <sup>1</sup>
Article Number	9010-0321	9010-0322
Performance Data Temperature		
Temperature range 10 °C above ambient temperature to [°C]	100	100
Temperature variation at 37 °C [± K]	0.3	0.3
Temperature fluctuation at 37 °C [± K]	0.1	0.1
Heating-up time to 37 °C [min]	15	15
Recovery time after 30 seconds door open at 37 °C [min]	4	4


<sup>1</sup> All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10%. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.



Description	BF720-230V <sup>1</sup>	BF720UL-240V <sup>1</sup>
Article Number	9010-0321	9010-0322
Electrical data		
Rated Voltage [V]	230	240
Power frequency [Hz]	50/60	60
Nominal power [kW]	1.75	1.85
Phase (Nominal voltage)	1~	1~
Measures		
Interior volume [L]	734	734
Net weight of the unit (empty) [kg]	165	165
Permitted load [kg]	315	315
Load per rack [kg]	45	45
Wall clearance back [mm]	160	160
Wall clearance sidewise [mm]	100	100
Internal Dimensions		
Width [mm]	1000	1000
Height [mm]	1300	1300
Depth [mm]	565	565
Doors		
Inner doors	2	2
Unit doors	2	2
Housing dimensions not incl. fittings and connections		
Width net [mm]	1165	1165
Height net [mm]	1590	1590
Depth net [mm]	870	870
Environment-specific data		
Energy consumption at 37 °C [Wh/h]	130	130
Sound-pressure level [dB(A)]	43	43
Fixtures		
Number of shelves (std./max.)	2/16	2/16

<sup>1</sup> All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10%. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

# DIMENSIONS incl. fittings and connections [mm]

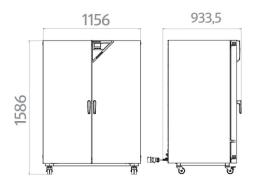



## **OPTIONS**

Designation	Description	*	ArtNo.
	left		
	10 mm	01	8012-1274
	30 mm	01	8012-1156
	50 mm	01	8012-1164
	right		
	10 mm	01	8012-1267
Access port with silicone plug	30 mm	01	8012-1153
pidg	50 mm	01	8012-1162
	top		
	10 mm	01	8012-126
	30 mm	01	8012-115
	50 mm	01	8012-115
	100 mm	01, 10	8012-116
Alarm function for overheating	Switchable acoustic alarm, with adjustable limit value on the independent temperature safety device	-	8012-163
Analog output 4-20 mA	for temperature values (output not adjustable)	02	8012-162
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test temperature	-	8012-111
	for temperature, measurement in center of chamber at specified temperature	-	8012-113
Calibration certificate,	temperature measurement incl. certificate and 27 measuring points at specified temperature	-	8012-159
temperature	temperature measurement incl. certificate, 15- 18 measuring points at specified temperature	-	8012-156
	temperature measurement incl. certificate, 9 measuring points at specified temperature	_	8012-154
Clock	battery-backed	-	8012-163
Ethernet interface	for Multi Management Software APT-COM™	-	8012-099
Inner chamber, reinforced	max. total load 350 kg	-	8012-182
Interior socket	Waterproof, switched flush-mounting box (CA3GD) for nominal voltage in unit interior, with cover and corresponding plug (nominal voltage; max. 500 W; max. 90 °C; protection class IP67)	07	8012-163
Pt 100 temperature sensor	additional flexible Pt 100, interior, for displaying the temperature on the unit display	-	8012-161
* Notes > See last page			

Notes > See last page

# **ACCESSORIES**

Designation	Description	*	ArtNo.
APT-COM™ 4 BASIC-	for simple logging and documentation requirements with up to 5 networked units.		
Edition	version 4, BASIC edition	19	9053-0039
APT-COM™ 4 GLP- Edition	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.		
Edition	version 4, GLP edition	19	9053-0042
APT-COM™ 4	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.		
PROFESSIONAL- Edition	version 4, PROFESSIONAL edition	19	9053-0040
Data Logger Kit	T 220: For continuous temperature logging from -90 °C to 220 °C. The kit includes 1 data logger, Pt 100 sensor with 2 m extension cable and 1 magnetic fixture for mounting to the BINDER unit	19	8012-0715
Data Logger Software	LOG ANALYZE software kit, configuration and evaluation software for all BINDER Data Logger Kits (incl. USB data cable)	19	8012-0821
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	-	1002-0016

<sup>\*</sup> Notes > See last page



Designation	Description	*	ArtNo.
	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature, CO <sub>2</sub> , O <sub>2</sub> , pressure, depending on unit		
	Digital in PDF format	-	7057-0001
O	Hard copy inside folder	-	7007-0001
Qualification documents	IQ/OQ/PQ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature, CO <sub>2</sub> , O <sub>2</sub> – or pressure, depending on unit		
	Digital in PDF format	-	7057-0005
	Hard copy inside folder	-	7007-0005
Rack	chrome plated	-	6004-0171
	stainless steel	-	6004-0155
Shelf, perforated	Stainless steel	_	6004-0187

<sup>\*</sup> Notes > See last page

## **SERVICES**

Designation	Description	*	ArtNo.
Installation services			
Unit installation	Unpacking and setting up of unit, connecting to existing connections, and checking function	13, 18	DL10-0100
Unit instructions	Unit function instructions for operation and programming of the controller	18	DL10-0500
Preventive maintenance			
Maintenance	Functional testing of all electrical and mechanical components, short calibration, documentation in the maintenance schedule	14, 18	DL20-0200
Calibration services			
	Expansion – including certificate, each additional measuring point in center of chamber at specified test temperature	14, 16, 17, 18	DL30-0102
Temperature calibration	Temperature calibration with 1 measuring point in center of chamber at a specified test temperature, including certificate	14, 16, 17, 18	DL30-0101
Temperature measurement, 18 measuring points	Temperature measurement with 18 measuring points at a specified test temperature, including certificate	14, 16, 17, 18	DL30-0118
Temperature measurement, 27 measuring points	Temperature measurement with 27 measuring points at a specified test temperature, including certificate	14, 16, 17, 18	DL30-0127
Temperature measurement, 9 measuring points	Temperature measurement with 9 measuring points at a specified test temperature, including certificate	14, 16, 17, 18	DL30-0109
Validation services			
Execution of IQ/OQ	Execution of IQ/OQ in accordance with qualification folder	15, 18	DL40-0100
Execution of IQ/OQ/PQ	Execution of IQ/OQ/PQ in accordance with qualification folder	15, 18	DL44-0500
Warranty service			
Warranty extension from 2 to 3 years	The warranty is extended from 2 to 3 years from the delivery date, wear parts are excluded	-	DL01-7041
Warranty extension from 2 to 5 years	The warranty is extended from 2 to 5 years from the delivery date, wear parts are excluded	-	DL01-7042

<sup>\*</sup> Notes > See last page



## **NOTES**

- Condensation may occur in the area around the access port. Access ports may be placed in custom locations for an additional charge.
- 02 UL mark is not granted when this option is used.
- Heat resistant only to max. 200 °C. 03
- Only available on units rated for 230 V.
- The additional heat input may affect the temperature behavior.
- Not in conjunction with the optional access port, door with window and interior lighting.
- 10 Not available on 23-liter units.
- Not available on 23- or 53-liter units.

  Only available on units rated for 230 V or 400 V. 11 12
- 13 Installation and connections take place at unit location; transport within the company only upon consultation.
- We recommend a BINDER service contract to cover unit inspections, calibrations and validations.
- OQ according to Yellow Paper = completed factory validation documentation of all OQ checklists.
- Sensor calibration is performed in an accredited calibration laboratory.
- Calibration is performed according to the BINDER factory standard. 17
- Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not 18 include a country-specific added fee (available on request).
- For additional accessories, refer to the Process documentation chapter.

BINDER GmbH Tuttlingen, Germany TEL +49 7462 2005 0 FAX +49 7462 2005 100 info@binder-world.com www.binder-world.com

BINDER Asia Pacific (Hong Kong) Ltd. Kowloon, Hong Kong, P.R. China TEL +852 39070500 FAX +852 39070507 asia@binder-world.com www.binder-world.com

**BINDER Environmental Testing** Equipment (Shanghai) Co., Ltd. Shanghai, P.R. China TEL +86 21 685 808 25 FAX +86 21 685 808 29 china@binder-world.com www.binder-world.com

Representative Office for CIS Moscow. Russia TEL +7 495 988 15 16 FAX +7 495 988 15 17 russia@binder-world.com www.binder-world.com

BINDER Inc. Bohemia, NY, USA TEL +1 631 224 4340 FAX +1 631 224 4354 usa@binder-world.com www.binder-world.us